



# SlimShaper PRO<sup>®</sup> BY ZARC

## DIRECTIONS FOR USE

SlimShaper PRO instruments for endodontic treatments:

- SlimShaper PRO Shaping Files (ZS1, ZS2, ZS3)

## COMPOSITION

The cutting surface of these instruments is made of a nickel-titanium alloy. DualWire<sup>®</sup> heat treated.

## 1) INDICATIONS FOR USE

These instruments have been designed to be used exclusively in a dental clinic or hospital setting by qualified users.

Application field: for the removal of dentin and shaping of the root canal.

## 2) CONTRAINDICATIONS

- As with all mechanically driven root canal instruments, SlimShaper PRO files should not be used in cases of severe and sudden apical curvatures due to heightened risk of separation.
- This product contains nickel and should not be used in patients with a known allergy to this material.

## 3) WARNINGS

- The files are single-use only, so they are not to be reused as this can cause deformation defects (bending, stretching), fractures, corrosion, loss of color or markings identification, etc., which would mean the instruments cannot meet the level of safety required for the intended use.
- Reuse increases the risk of cross-contamination and breakage.
- There is an increased risk of file separation when used after several disinfection or sterilization cycles.

## 4) PRECAUTIONS FOR USE

- Straight-line access is a prerequisite for proper root canal treatment, and with SlimShaper PRO files there is no exception.
- These instruments should not be immersed in a sodium hypochlorite solution.

- Clean files frequently during instrumentation, inspecting for any signs of distortion or wear, such as missing flutes or nicks.
- Irrigate frequently, recapitulate and irrigate the canal throughout the whole procedure.
- SlimShaper PRO files should only be used in regions of the canal that have a confirmed and reproducible glide path. Establish a reproducible glide path using hand files of at least an ISO 015 size.
- SlimShaper PRO files are manufactured with a process that results in a gold, blue or pink appearance. Due to this process, SlimShaper PRO files may appear slightly curved. This is not a manufacturing defect. While the file can be easily straightened using only your fingers, it is not necessary to straighten the file prior to use. Once inside the canal, SlimShaper PRO files will follow the anatomy.
- Always use minimal apical pressure. Never force the files inside the canal.
- This product should be disposed of as medical waste.
- For optimal usage, torque control devices are recommended.
- SlimShaper PRO rotary files can be used at motor speeds between 300 rpm and 500 rpm.

Recommended motor settings:

SlimShaper PRO		
File size	Speed (rpm)	Torque (Ncm)
SlimShaper PRO ZS1, ZS2 and ZS3	500	3

The speed and torque settings indicated in the table above are demonstrative and can be modified according to user preferences and equipment possibilities.

## 5) ADVERSE REACTIONS

SlimShaper PRO is not recommended for use in patients with a known allergy to nickel. The use of this product in these patients may cause: difficulty breathing, swelling of the face or eyes, hives or rash. If any of these symptoms occur, patients should be advised to contact a dental professional immediately.

## 6) STEP-BY-STEP INSTRUCTIONS FOR USE OF SlimShaper PRO FILES

### 6.1 Radiographic assessment

Review several horizontal angle radiographs to diagnostically determine the width, length, and curvature of the canals.

### 6.2 Access preparation

Create a straight-line access to the canal entrance, keeping in mind the taper, flattening, and finishing of the internal axial walls.

### 6.3 SlimShaper PRO SHAPING TECHNIQUE



The crown down technique is the technique of choice for rotary instruments:

- Create straight-line access to the canal entrance orifice.
- In the presence of NaOCl, perform an exploration of the coronal 2/3 with hand files 10 and 15. Gradually work these instruments until a reproducible glide path is confirmed. Alternatively, mechanized glide path files (such as Z-Glider®) can be used after a 10 hand file.
- In the presence of NaOCl, “float” the ZS1 in the canal and passively “follow” the glide path. Before light resistance is encountered, “brush” laterally cutting the dentin on the outstroke to improve straight-line access and apical progression. Always brush away from the furcation.
- Continue shaping with ZS1 as described until the depth of a 15 hand file is reached.
- Use ZS2, exactly as described for ZS1, until the depth of a 15 hand file is reached.
- In the presence of a viscous chelator or NaOCl, explore the apical third with hand files 10 and 15 and gradually work them until they are loose at length.
- Establish working length, confirm patency and verify the presence of a smooth and reproducible glide path in the apical third.
- Do not use brushing action.
- Reconfirm the working length, irrigate, recapitulate and reirrigate, especially in more curved canals.
- Use ZS3 finishing file without brushing action, going deeper with each insertion than the previous insertion until working length is reached, spending the least amount of time possible at this length. Do not leave the file at working length for longer than one second. Reach it and withdraw it.
- Gauge the foramen with a 20 hand file. If the instrument has a close fit at length, the canal is properly shaped and ready for obturation.
- The SlimShaper PRO sequence is the same regardless of the length, diameter or curvature of the canal.

## 7) CLEANING, DISINFECTION AND STERILIZATION

Dental instrument reprocessing procedure.

### I - FOREWORD

**Instruments marked “sterile” do not require any specific treatment before first use.**

**For the rest of the instruments that are not marked “sterile”, cleaning and sterilization is necessary before using them for the first time in accordance with section III – STEP-BY-STEP INSTRUCTIONS of these directions for use.**

### Excluded devices:

Uniclip and Mooser calcinable plastic posts cannot be sterilized and must be disinfected by immersion in NaOCl (2.5 % at least) during 5 minutes at ambient temperature.

### II - GENERAL RECOMMENDATIONS

- 1) Use only a detergent solution with a disinfectant action approved for its effectiveness (VAH/DGHM-listing, CE marking, FDA approval) and follow the instructions of the manufacturer of said product. Anticorrosive cleaning products and disinfectants are recommended for all metal instruments.
- 2) For your own safety, please use personal protective equipment (gloves, goggles and mask).
- 3) The user is responsible for the cleaning and sterilization of the product for the first



cycle and each subsequent use, as well as the use of dirty or damaged instruments after sterilization.

- 4) The safest option for the professional is to use our instruments only once. After each use, they must be thoroughly inspected before being used again; the presence of defects such as deformations (bent, stretched), fractures, corrosion, loss of color identification, marking, etc. are indications that the instruments cannot meet the level of safety required for their intended use and should therefore be discarded.
- 5) Instruments marked as single-use only are not approved for reuse.
- 6) For the final rinse it is mandatory to use deionized water, either using an automatic washer-disinfector or a manual cleaning method. For the other rinsing steps tap water use is allowed.
- 7) Instruments with plastic handles, and NiTi instruments should not be used with Hydrogen Peroxide (H<sub>2</sub>O<sub>2</sub>) solution which is known to degrade them.
- 8) Only the active part of the NiTi instrument that is in contact with the patient should be immersed in a NaOCl solution with a concentration NOT greater than 5%.
- 9) Prevent the instrument from drying out before or during pre-disinfection or cleaning. Dried biological material can be difficult to remove.
- 10) Do not place identifying labels or markers directly on the instrument.

Operation		Description and Warnings
1	<b>Disassembling</b>	<ul style="list-style-type: none"> <li>• Remove and discard the silicone stoppers.</li> </ul>
2	<b>Rinsing</b>	<ul style="list-style-type: none"> <li>• Rinse abundantly (for at least 1 minute) under running water at ambient temperature. While rinsing, use a soft brush (nylon, polypropylene, acrylic) for pre-cleaning until visible impurities are removed.</li> </ul>
3a	<b>Automated cleaning with washer-disinfector</b>	<ul style="list-style-type: none"> <li>• Place the instrument in a kit, support, or container made of stainless steel or titanium.</li> <li>• Execute the defined cycle with detergent solution (for example, Metrex EmPowder concentration 1:128 - 1:512) for at least 5 minutes in the washer-disinfector at 20°C ~ 40°C).</li> </ul>
3b.i	<b>Manual cleaning assisted by an ultrasonic device</b>	<ul style="list-style-type: none"> <li>• Place the instrument in a kit, support, or container made of stainless steel, polypropylene, or titanium.</li> <li>• Submerge it in a detergent solution (for example, Metrex EmPowder concentration 1:128) with cleaning properties. If applicable, soak it for at least 15 minutes with the help of an ultrasonic device.</li> <li>• Rinsing: perform a long rinse (at least 1 minute) under running deionized water at a temperature of 20°C ~ 40°C.</li> <li>• Drying: dry with a single-use non-woven cloth or with a hot air dryer no hotter than 110°C.</li> </ul>



Operation		Description and Warnings
3b.II	<b>Manual disinfection with washer-disinfector</b>	<ul style="list-style-type: none"> <li>Place the instrument in a kit, support, or container made of stainless steel or titanium.</li> <li>Execute the defined cycle with a mild neutral enzymatic cleaning agent solution (for example, Metrex EmPowder concentration 1:128) for at least 5 minutes in the washer-disinfector at a temperature of &gt;90°C, AO &gt;3000.</li> <li><b>Note:</b> <ol style="list-style-type: none"> <li>Discard instruments with obvious defects (broken, bent, etc.).</li> <li>When the instruments are placed in a cleaning kit, support, or container, avoid any contact between them.</li> <li>Follow instructions and concentrations given by the manufacturer of the detergent solution (see also general recommendations).</li> <li>Follow the instructions of the washer-disinfector and verify the success criteria after each cycle have been met as stated by the manufacture.</li> <li>The final rinse step should be with deionized water. For other steps follow the water quality defined by the manufacturer. Place the instrument in a kit, support, or container made of stainless steel or titanium to avoid any contact between devices or posts.</li> </ol> </li> </ul>
3b.III	<b>Rinsing</b>	<ul style="list-style-type: none"> <li>Rinse abundantly (for at least 1 minute) under running water at ambient temperature.</li> <li>Use deionized water for rinsing.</li> <li>If the previously used cleaning solution contains a corrosion inhibitor, rinsing is recommended just before autoclaving.</li> </ul>
3b.IV	<b>Drying</b>	<p>Devices should be carefully dried before inspection and packaging</p> <ul style="list-style-type: none"> <li>Dry with a single-use non-woven cloth or with a hot air dryer at a maximum temperature of 110 °C.</li> <li>The devices should be dried until visual traces of moisture are eliminated.</li> <li>Particular attention has to be paid to effectively dry joints or cavities inside the device.</li> </ul>
4	<b>Inspection</b>	<ul style="list-style-type: none"> <li>Inspect the operation of the devices.</li> <li>Inspect devices and identify those with defects. <ol style="list-style-type: none"> <li>Dirty devices need to be cleaned again.</li> <li>Do not reuse the silicone stoppers.</li> <li>Dispose of any defective devices.</li> </ol> </li> </ul>



	Operation	Description and Warnings
5	<b>Packaging</b>	<p>Place the devices in a kit, support or container to avoid any contact between instruments or posts and pack the devices in “Sterilization Pouches”. (Use packaging that is resistant to a temperature of 141 °C (286 °F) and that is UNE-EN ISO 11607 compliant).</p> <ol style="list-style-type: none"> <li><b>1)</b> Avoid any contact between instruments or posts during sterilization. Use kits, supports or containers.</li> <li><b>2)</b> For sharp devices not contained within a box, silicone tubing should be placed around the devices to prevent puncture of the packaging.</li> <li><b>3)</b> Seal pouches per the pouch manufacturer’s recommendation. If a thermo-sealer is used, the process must be validated.</li> <li><b>4)</b> Check the expiration date of the pouch indicated by the pouch manufacturer to determine the shelf life of the sterile product.</li> </ol>
6	<b>Sterilization</b>	<ul style="list-style-type: none"> <li>• For these devices, steam sterilization is recommended at 132°C / 273°F for 4 minutes, for the purpose of de-activating potential prions.</li> <li>• Instruments and posts must be sterilized according to the package labelling.</li> <li>• Place pouches in steam sterilizer per sterilizer manufacturer’s recommendation.</li> <li>• Use only steam sterilizers that meet the requirements of EN 13060 (class B, small sterilizer), EN 285 (full size sterilizer).</li> <li>• Use a validated sterilization procedure according to ISO 17665 with a minimum drying time of 20 minutes.</li> <li>• Follow the sterilizer maintenance procedure indicated by the sterilizer manufacturer.</li> <li>• Control the efficiency and the acceptance criteria of the sterilization procedure (packaging integrity, no humidity, no color change of the packaging, positive physicochemical indicators, conformity of actual cycle parameters to reference cycle parameters).</li> <li>• Store traceability records and define shelf life according to packaging manufacturer guidelines.</li> <li>• Shorter sterilization cycles according to local regulations are possible but are not guaranteed to de-activate prions.</li> </ul>
7	<b>Storage</b>	<p>Keep devices in sterilized packaging in a clean environment, away from sources of moisture and direct sunlight. Store at ambient temperature.</p> <ol style="list-style-type: none"> <li><b>1)</b> Sterility cannot be guaranteed if the packaging is open, damaged or wet.</li> <li><b>2)</b> Check the packaging and medical devices before use (integrity of the packaging, absence of humidity and expiration date).</li> </ol>



## 8) TRANSPORT

- To prevent the medical device from being damaged during transport, the use of specific racks, trays or rigid containers may be recommended.
- When the package is broken, it is no longer sterile. It must be cleaned, disinfected and sterilized before use.

## 9) STORAGE AND EXPIRATION

- Avoid storing in places with high temperatures, humidity and/or direct sunlight. Keep liquids away. Store at ambient temperature.
- Do not damage or puncture the packaging materials.
- This product is subject to improvement without prior notice. Apply the “first in, first out” approach to inventory management.
- Do not store under a germicidal lamp to prevent deterioration.
- The shelf life of endodontic files is 5 years.




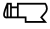






## 10) DISPOSAL

- This product should be disposed of as medical waste.
- For correct disposal, always respect national laws and recommendations from authorities.






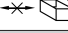


## 11) PACKAGING


- Packing configurations: 3pcs/pack or 6pcs/pack in aluminium foil box.
- Assortment: 1 piece or 2 pieces of each size in a single pack.

## 12) IDENTIFICATION OF RELATED SYMBOLS

Symbols	Identification
 FOR DENTAL USE ONLY	Product designated for dental use only
	Do not reuse – Single-use only
	Workpiece Material: Nickel Titanium
	Rotating handle
	Medical Device
Heat activation	Pre-bendable
	Recommended rotation speed
	Recommended torque for use
	Batch number
	Sterilized by radiation
	Reference number



Symbols	Identification
	Expiration date
	Consult directions for use
	Authorized representative in the EU
	Manufacturer
	Autoclavable at the specified temperature
	Non-returnable if the seal is broken
	Stopper material: silicone
	Do not use if packaging is damaged. Consult directions for use.

 **Mfg. for NuSmile, Ltd.**  
 3315 West 12th St.  
 Houston, TX 77008 USA  
 nusmile.com - 713.861.0033

All rights reserved.

**SlimShaper PRO<sup>®</sup>**

FOR DENTAL USE ONLY  
**IFU-130 Rev 0**

